

Product information for suction handpieces

Product description

Dürr Dental offers various suction handpieces for the large suction hose (spray-mist suctioning) and the small suction hose (saliva ejector).

- Large suction handpiece with diameter 11 and 16 mm
- Small suction handpiece in angled or straight form
- With slide and/or rotating bar to regulate suction performance
- Ball joint for ergonomic work

Material

PA

Slide / rotating bar: TPU

O-ring: silicone

Colour

RAL 7042

Material durability

Plastics and O-rings are durable against:


- Orotol plus and MD555
- FD 350 Disinfecting tissues for preliminary cleaning
- ID 215 Enzymatic instrument cleaner
- ID 213 Instrument disinfection

All handpieces can be cleaned in the RDG and steam-sterilised in an autoclave at 134°C (for information on manual and machine processing, see instructions, order no. 9000-606-18)


Suction of spray mist

Ø 16 mm


Ø 11 mm

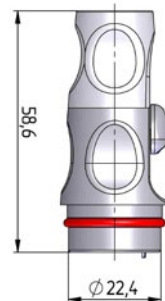


7600A010-00

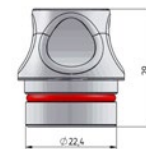

7600A040-00

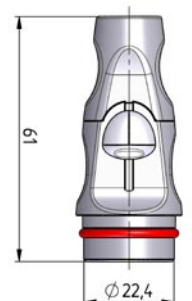


7600A010-75

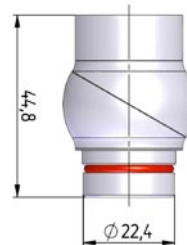


7600A030-00


Saliva ejector

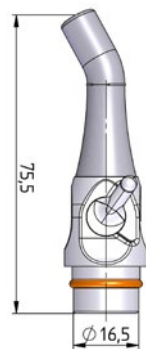


7600A020-00

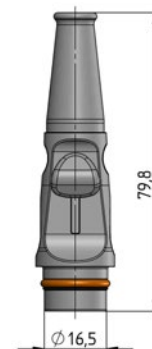

7600A025-00

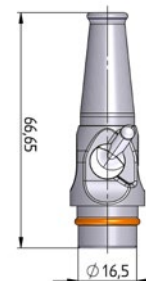

7600A020-55

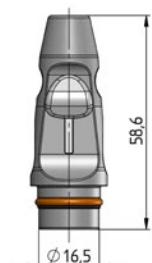

7600A025-50